



CORONAVIRUS TESTING EXPLAINED



TEST TYPES

ANTIBODY TESTING

Use Cases

- Serological Testing and Instant Screening

Result wait time

- Instant: 10 Minutes

Sample Collection

- Serum, Plasma or Whole Blood

Method of detection

- Presence of Antibodies in Blood Sample

www.statnews.com

www.sciencemag.org



PCR TESTING

Use Cases

- Gold Standard for Point of Care Diagnostic

Result Wait Time

- Instant: 5-45 Minutes
- Lab: 1-10 Days

Sample Collection

- Nasal Swab

Method of Detection

- Partial Genetic Markers of the Virus

Downsides

- Requires Expensive PCR Machine

www.statnews.com

www.cnn.com


www.wired.com

www.time.com



CORONAVIRUS ANTIBODIES

COVID-19 STAGES, TARGET TIMELINES AND CLINICAL

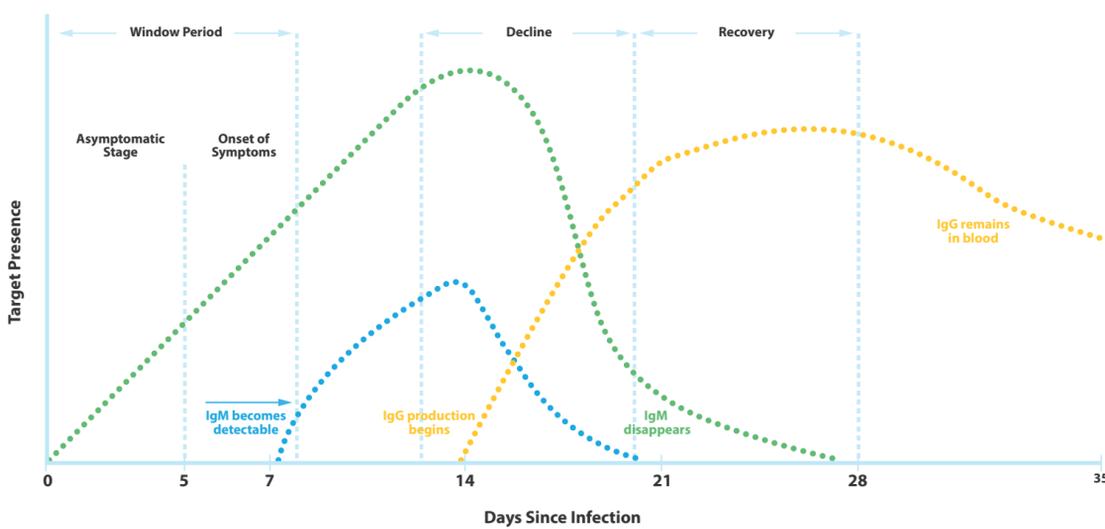

- SARS-CoV-2 RNA
- IgM Antibody
- IgG Antibody

DISCLAIMER:

This is for information purposes only. Detection windows will vary from person to person.

Test Results			Clinical Significance
RNA	IgM	IgG	
+	-	-	Patient may be in the window period of the infection
+	+	-	Patient may be in the early stage of infection
+	+	+	Patient is in the active phase of infection
+	-	+	Patient may be in the late or recurrent stage of infection
-	-	+	Patient may have had a past infection, and has recovered



DIFFERENCE BETWEEN IGG AND IGM

IgM Only

- Patient Is Fighting Infection

Immunoglobulin M (IgM), which is found mainly in the blood and lymph fluid, is the first antibody to be made by the body to fight a new infection.

IgG Only

- Patient Has Recovered

Immunoglobulin G (IgG), the most abundant type of antibody, is found in all body fluids and protects against bacterial and viral infections.

Both

- Patient is still fighting infection

kidshealth.org



COVID-19 "CORONAVIRUS" IGG/IGM RAPID TEST KIT

Detection Window

- IgM may appear in blood within 3-5 days. IgG will appear as soon as 1-2 weeks

More information, such as the product insert, clinical trial results, and CE documentation is available at request.

Sample

- Test can work with whole blood, plasma and serum samples

To request more information about our Rapid Coronavirus (Covid-19) test kit-such visit:

Storage

- The kit can be stored at room temperature or refrigerated (2-30°C)

Shelf Life

- 24 months from manufacture date



My CCM Safe Certified PPE
Moe@mychroniccaremanagement.net
www.mychroniccaremanagement.net

COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma)

INTENDED USE

COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma) is a solid phase immunochromatographic assay for the rapid, qualitative and differential detection of IgG and IgM antibodies to 2019 Novel Coronavirus in human whole blood, serum or plasma. This test provides only a preliminary test result. Therefore, any reactive specimen with the COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma) must be confirmed with alternative testing method(s) and clinical findings.

INTRODUCTION

Coronaviruses are enveloped RNA viruses that are distributed broadly among humans, other mammals, and birds and that cause respiratory, enteric, hepatic, and neurologic diseases. Seven coronavirus species are known to cause human disease. Four viruses - 229E, OC43, NL63, and HKU1 - are prevalent and typically cause common cold symptoms in immunocompetent individuals. The three other strains - severe acute respiratory syndrome coronavirus (SARS-CoV), Middle East respiratory syndrome coronavirus (MERS-CoV) and 2019 Novel Coronavirus (COVID-19) - are zoonotic in origin and have been linked to sometimes fatal illness. IgG and IgM antibodies to 2019 Novel Coronavirus can be detected with 2-3 weeks after exposure. IgG remains positive, but the antibody level drops overtime.

PRINCIPLE

The COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma) is a lateral flow immunochromatographic assay. The test uses anti-human IgM antibody (test line IgM), anti-human IgG (test line IgG) and goat anti-rabbit IgG (control line C) immobilised on a nitrocellulose strip. The burgundy colored conjugate pad contains colloidal gold conjugated to recombinant COVID-19 antigens conjugated with colloid gold (COVID-19 conjugates) and rabbit IgG-gold conjugates. When a specimen followed by assay buffer is added to the sample well, IgM &/or IgG antibodies if present, will bind to COVID-19 conjugates making antigen antibodies complex. This complex migrates through nitrocellulose membrane by capillary action. When the complex meets the line of the corresponding immobilized antibody (anti-human IgM &/or anti-human IgG) the complex is trapped forming a burgundy colored band which confirm a reactive test result. Absence of a colored band in the test region indicates a non-reactive test result.

The test contains an internal control (C band) which should exhibit a burgundy colored band of the immunocomplex goat anti rabbit IgG/rabbit IgG-gold conjugate regardless of the color development on any of the test bands. Otherwise, the test result is invalid and the specimen must be retested with another device.

MATERIALS SUPPLIED

25 sealed pouches each containing a test cassette, a dropper and a desiccant
1 Buffer
1 package insert

MATERIAL REQUIRED BUT NOT PROVIDED

1. Specimen collection containers
2. Lancets (for fingerstick whole blood only)
3. Centrifuge (for plasma only)
4. Timer
5. Heparinized capillary tubes and dispensing bulb (for fingerstick whole blood only)

STORAGE AND STABILITY

The kit can be stored at room temperature or refrigerated (2-30°C). The test device is stable through the expiration date printed on the sealed pouch. The test device must remain in the sealed pouch until use. DO NOT FREEZE. Do not use beyond the expiration date.

WARNINGS AND PRECAUTIONS

1. For professional In Vitro diagnostic use only. Do not use after expiration date.
2. This package insert must be read completely before performing the test. Failure to follow the insert gives inaccurate test results.
3. Do not use it if the tube/pouch is damaged or broken.
4. Test is for single use only. Do not re-use under any circumstances.
5. Handle all specimens as if they contain infectious agents. Observe established precautions against microbiological hazards throughout testing and follow the standard procedures for proper disposal of specimens.
6. Wear protective clothing such as laboratory coats, disposable gloves and eye protection when specimens are assayed.
7. Humidity and temperature can adversely affect results.
8. Do not perform the test in a room with strong air flow, ie. electric fan or strong air-conditioning.

SPECIMEN COLLECTION

1. COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma) can be performed using either whole blood, serum or plasma.
2. Separate serum or plasma from blood as soon as possible to avoid hemolysis. Use only clear, non-hemolyzed specimens.
3. Testing should be performed immediately after specimen collection. Do not leave the specimens at room temperature

for prolonged periods. Serum and plasma specimens may be stored at 2-8°C for up to 3 days. For long term storage, specimens should be kept below -20°C. Whole blood collected by venipuncture should be stored at 2-8°C if the test is to be run within 2 days of collection. Do not freeze whole blood specimens. Whole blood collected by fingerstick should be tested immediately.

4. Bring specimens to room temperature prior to testing. Frozen specimens must be completely thawed and mixed well prior to testing. Specimens should not be frozen and thawed repeatedly.
5. If specimens are to be shipped, they should be packed in compliance with local regulations covering the transportation of etiologic agents.

TEST PROCEDURE

Allow test cassette, specimen, buffer and/or controls to equilibrate to room temperature (15-30°C) prior to testing.

1. Remove the test cassette from the sealed foil pouch and use it as soon as possible. Best results will be obtained if the assay is performed within one hour.
2. Place the test device on a clean and level surface.

For Serum or Plasma Specimens:


With a 5 µL mini plastic dropper provided, draw serum/plasma specimen to exceed the specimen line as showed in the following image and then transfer drawn serum/plasma specimen into the sample well (S). Then add 2 drops (about 80 µL) of sample buffer to the buffer well (B) immediately. Avoid air bubbles.

Note: Practice a few times prior to testing if you are not familiar with the mini dropper. For better precision, transfer specimen by pipette capable to deliver 5 µL of volume.

For Whole Blood Specimen:

Hold the 5 µL mini plastic dropper vertically and transfer 1 drop of whole blood (about 10 µL) to the specimen well(S) of the test device, then add 2 drops (about 80 µL) of sample buffer to the buffer well (B) immediately. Avoid air bubbles.

3. Wait for the colored line(s) to appear. The result should be read in 10 minutes. Positive results may be visible as soon as 2 minutes. Do not interpret the result after 15 minutes.

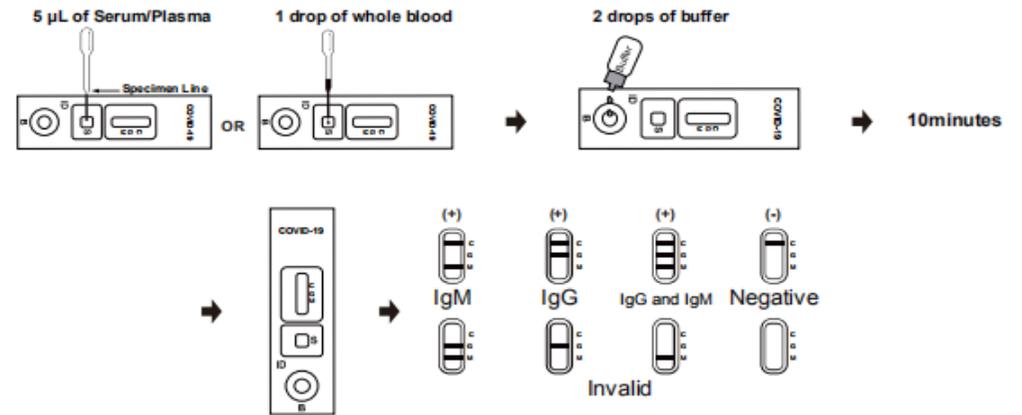

INTERPRETATION OF RESULTS

NEGATIVE: If only the C band is present, the absence of any burgundy color in the both T bands (IgG and IgM) indicates that no anti-COVID-19 antibodies are detected in the specimen. The result is negative.

IgM POSITIVE:

In addition to the presence of C band, if only IgM band is developed, the test indicates for the presence of IgM anti-COVID-19 in the specimen. The result is IgM anti-COVID-19 positive.

IgG POSITIVE:

In addition to the presence of C band, if only IgG band is developed, the test indicates for the presence of IgG anti-COVID-19 in the specimen. The result is IgG anti-COVID-19 positive.

IgG and IgM POSITIVE:

In addition to the presence of C band, both IgG and IgM bands are developed, the test indicates for the presence of both

IgG and IgM anti-COVID-19 in the specimen. The result is IgG and IgM anti-COVID-19 positive.

INVALID:

Control line fails to appear. Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for control line failure. Review the procedure and repeat the test with a new test cassette. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

QUALITY CONTROL

A procedural control is included in the test. A red line appearing in the control region (C) is the internal procedural control. It confirms sufficient specimen volume and correct procedural technique.

Control standards are not supplied with this kit; however it is recommended that positive and negative controls be tested as a good laboratory practice to confirm the test procedure and to verify proper test performance.

LIMITATIONS

3. Use fresh samples whenever possible. Frozen and thawed samples (especially repeatedly) contain particles that can block the membrane. This slows the flow of reagents and can lead to high background color, making the interpretation of results difficult.

4. Optimal assay performance requires strict adherence to the assay procedure described in this insert sheet. Deviations may lead to aberrant results.

5. A negative result for an individual subject indicates absence of detectable anti-COVID-19 antibodies. However, a negative test result does not preclude the possibility of exposure to or infection with COVID-19.

6. A negative result can occur if the quantity of the anti-COVID-19 antibodies present in the specimen is below the detection limits of the assay, or the antibodies that are detected are not present during the stage of disease in which a sample is collected.

7. Some specimens containing unusually high titer of heterophile antibodies or rheumatoid factor may affect expected results.

8. As with all diagnostic tests, a definitive clinical diagnosis should not be based on the result of a single test, but should only be made by the physician after all clinical and laboratory findings have been evaluated.

PERFORMANCE CHARACTERISTICS

1. Clinical Performance for IgM Test

The samples from susceptible subjects were tested by The COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma) and by a commercial IgM EIA kit.

Relative Sensitivity: 95.7%, Relative Specificity: 97.3%, Overall Agreement: 96.8%

2. Clinical Performance for IgG Test

The samples from susceptible subjects were tested by the COVID-19 IgG/IgM Rapid Test Cassette (Whole Blood/Serum/Plasma) and by a commercial IgG EIA kit.

Relative Sensitivity: 91.8%, Relative Specificity: 96.4%, Overall Agreement: 95.0%

REFERENCE

1. Weiss SR, Leibowitz JL. Coronavirus pathogenesis. *Adv Virus Res* 2011; 81: 85-164.
2. Masters PS, Perlman S. Coronaviridae. In: Knipe DM, Howley PM, eds. *Fields virology*. 6th ed. Lippincott Williams & Wilkins, 2013: 825-58.
3. Su S, Wong G, Shi W, et al. Epidemiology, genetic recombination, and pathogenesis of coronaviruses. *Trends Microbiol* 2016; 24: 490-502.
4. Cui J, Li F, Shi ZL. Origin and evolution of pathogenic coronaviruses. *Nat Rev Microbiol* 2019; 17: 181-192.